



## 50/15 PRIME M

**Cod. PR-OF200-040-M**

Oil free reciprocating air compressors, 2 cylinders, with a compressed air FAD of 152 l/min @ 5 bar pressure. Receiver capacity 40 L, with internal surface treatment to protect against corrosion. Complete with membrane dryer to ensure maximum air purity. This compressor is available also in **10 BAR VERSION**, visit the appropriate section!


 Maximum simultaneous use: 2

 Dental Chairs: 3

Medical Device CE 0425

Cod.	Air Dryer	Air capacity	FAD Fad @ 5bar	Motor	Voltage Hz	Receiver	Noise Level	RPM	Duty Cycle	max P.	Size	net Weight
PR-OF200-040-M	Yes	m <sup>3</sup> /min/l/min 0.24 240	m <sup>3</sup> /min l/min 0.152 152	kW 1.5	230V 50Hz	l 40	db(A) 66	1360	70 %	bar 7	mm 710x410x770	Kg 50
PR-OF200-040-M-400	Yes	m <sup>3</sup> /min/l/min 0.255 255	m <sup>3</sup> /min l/min 0.152 152	kW 1.5	400V 50Hz	l 40	db(A) 66	1470	70 %	bar 7	mm 710x410x770	Kg 50
PR-OF200-040-M-6-UL	Yes	m <sup>3</sup> /min/l/min 0.28 280	m <sup>3</sup> /min l/min 0.165 165	kW 1.5	220V 60Hz	l 40	db(A) 66	1760	70 %	bar 7	mm 710x410x770	Kg 50
PR-OF200-040-M-110	Yes	m <sup>3</sup> /min/l/min 0.28 280	m <sup>3</sup> /min l/min 0.165 165	kW 1.5	110V 60Hz	l 40	db(A) 66	1760	70 %	bar 7	mm 710x410x770	Kg 50

## Air delivery performance

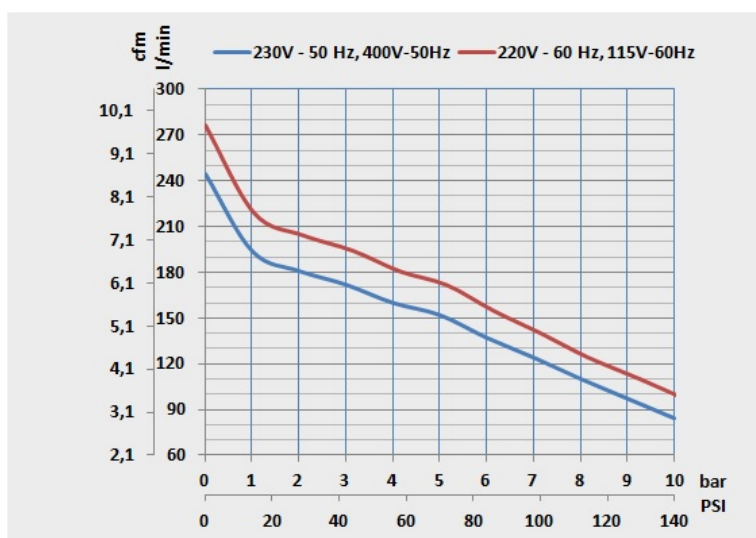

## Maintenance kits

Kit.	Description	Type	Hours
------	-------------	------	-------

TA-FL-026-K-SCP	<p><b>PRIME 2-3 M filter cartridges kit</b></p> <p>Cartridge replacement kit for PURE AIR PRIME 2,3 M HP range (cartridges + suction filter + membrane filter dryer cartridges). Replace the filtration cartridges every 2 years.</p>	Ordinary	1.500h/2y
PR-GC-031-K-SCP	<p><b>PRIME 1,5 kW</b></p> <p>The PRIME 1.5 kW spare parts kit includes: intake filters, piston rings, piston ring, O-ring, line filter cartridge, piston guide band and complete valve plate. The recommended maintenance intervals for this compressor are 4000 hours / 4 years.</p>	Extraordinary	4.000h/4y